


Adult Sternal Retractor

**Unique versatile
system minimizes
patient complications
and maximizes valve
exposure**



- Versatile System
- Innovative Design
- Minimizes Patient Complications

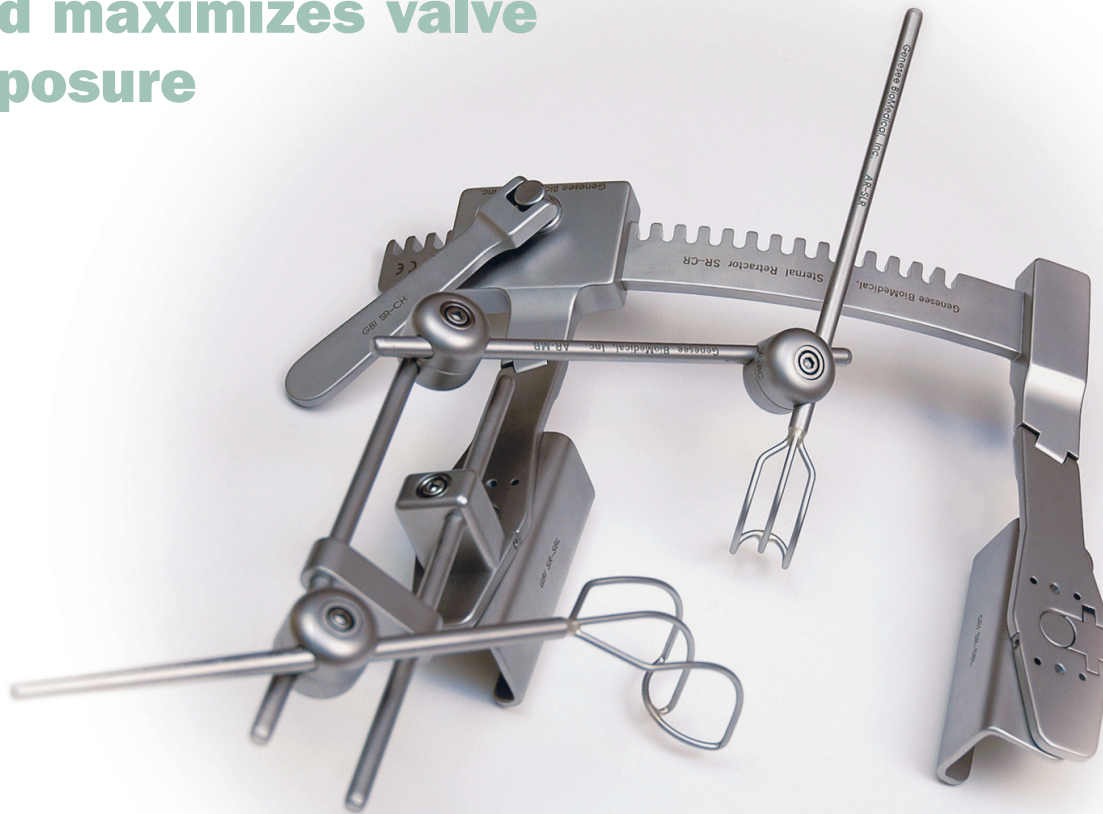

Genesee BioMedical, Inc.

Peak Performance Products for Cardiac Surgery

1.800.786.4890

www.geneseebiomedical.com



Adult Sternal Retractor with Companion Atrial & Aortic Retractors

A Multipurpose Retractor for Valve Replacement or Repair and CABG On or Off Pump

Versatile System


- Platform for Off-Pump Beating Heart Surgery
- Sternal Retractor for On-Pump Coronary Artery Bypass Grafting (CABG)
- Valve Retractor Maximizes Exposure for Valve Replacement or Repair

Innovative Design

- Curved Rack with Hinged Arms Allows Retractor to Hug Patient's Body Unobtrusively
- Unique Ball Clamps, Rotary H Bar and Pillar Clamp Provide Fast Setup and Unlimited Versatility in Tissue Retraction
- Changeable Blades Allow Surgeon to Adjust Retractor to Patient Anatomy
- Unique Shape of Ball and Pillar Clamps Prevent Suture Snagging
- Interchangeable Rakes Combined with Ball Clamps and H Bar Provide Excellent Exposure for Valve Replacement and Repair
- Simple Latch System Allows Surgeon to Rapidly Change Pivoting Blades to Patient Size (optional)

Minimizes Patient Complications

- Detachable Pivoting Retractor Blades Align with the Edges of the Sternum, Applying Pressure Evenly Across the Sternum
- Unique Blade Design Minimizes First Rib and Sternal Fractures
- Curved Rack Allows Sternum to be Opened in Accordance with the Natural Curvature of the Rib Cage, Minimizing Sternal Damage
- Curved Blades (optional)

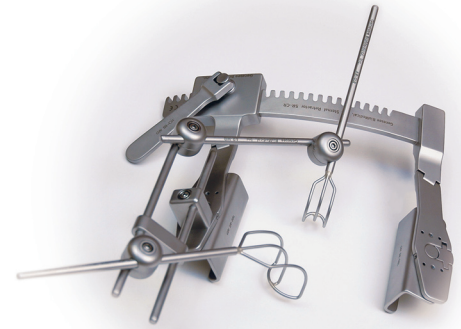

Ordering Information

Catalog Number	Product Description
SRC-AR	Sternal Retractor (Curved Rack) with Atrial & Aortic Retractors
SRC-1	Sternal Retractor Only
ST-ST5	Standard Sterilization Tray
SRCL-AR	Sternal Retractor with Latch System and Atrial & Aortic Retractors
SRCL-1	Sternal Retractor with Latch System Only

